

# High image quality in a compact laser imager.



## Optimal diagnostic images

With DRYPRO Sigma II, create images with unrivaled clarity and sharpness. Konica Minolta brings to you the latest precision optics which produce a 50 µm pixel pitch, for general radiography and mammography needs. Powerful image-processing algorithms are used to simultaneously optimize image smoothness and text sharpness. Diagnostic clarity is preserved and patient data is always legible, regardless of the print size.




Konica Minolta DRYPRO Sigma II

- Stable finished quality with automatic measurement to maintain the output density using the density control function
- User friendly design allows for easy operation and intuitive usability
- Powerful image-processing algorithms optimize both the image smoothness and the text sharpness
- Compact footprint dry imager with 2 trays and a processing capacity of 70 sheets/hr. in 14"x17" film

# Flexible imaging that meets your needs, in any environment.

## Space saving design

DRYPRO Sigma II is a desktop 2-ch dry imager. It features a compact footprint as small as 65 x 63 cm<sup>2</sup> and is designed to integrate into various hospital and clinic environments.

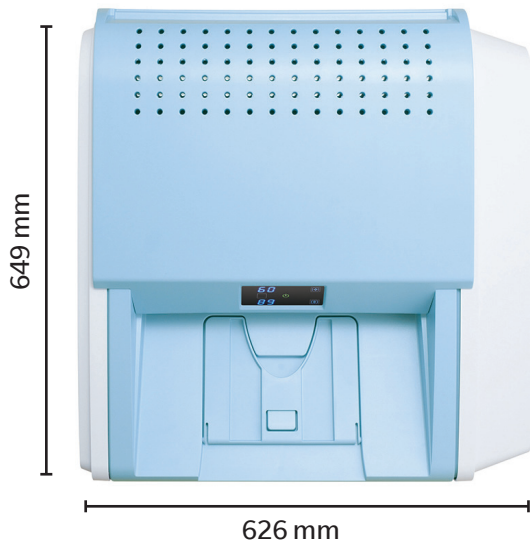

## Flexibility in a film daylight settings

SD-S Film is available in 4 sizes and two film trays. The two film trays can be used simultaneously making the system highly versatile. The system also calibrates whenever a tray of film with a new lot-number is loaded.



## DRYPRO Sigma II package


The DRYPRO Sigma II package is easy to purchase and use. Konica Minolta provides the best solutions for hospitals and clinics, and offers a variety of products to suit your needs.



## Network functions designed for open and flexible environments

DRYPRO Sigma II is a networked imager capable of connecting directly to DICOM print compliant devices.

In addition to DICOM basic grayscale print functions, DRYPRO Sigma II supports Presentation LUT, which enables printed film to more accurately match diagnostic monitors.

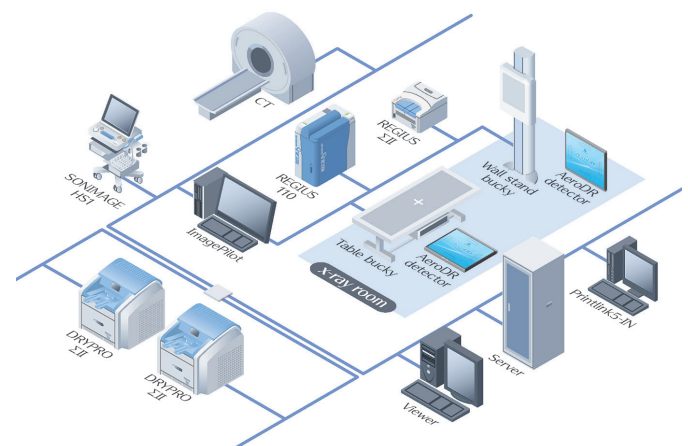

DRYPRO and AeroDR are registered trademarks of Konica Minolta Healthcare Americas, Inc.  
© 2021 Konica Minolta Healthcare Americas, Inc.

Sample workflow set-up



**Konica Minolta Healthcare Americas, Inc.**  
411 Newark Pompton Turnpike  
Wayne, New Jersey 07470  
Tel: +1 (973) 633-1500  
km.marketing@konicaminolta.com  
[healthcare.konicaminolta.us](http://healthcare.konicaminolta.us)

M1788 0121 RevA